

Discovery Hybrid Rheometer



Site Preparation Guide

Table of Contents


Table of Contents	2
Ideal Setup	3
System Components.....	4
Instrument Measurements	5
Utility Requirements.....	6
Computer Requirements	7-8
Hardware	7
Software.....	8
Temperature Systems.....	9-11
EHP, ETC, Peltier Plate/Peltier Concentric Cylinder	9
UHP	10-11
Site Preparation Checklist	12
TA Instrument Offices.....	13

Ideal Setup



IDEAL PLACEMENT AND BENCH MEASUREMENTS

Select a location with adequate floor space and a rigid laboratory bench that is level and is in a vibration-free environment. A marble table is recommended.

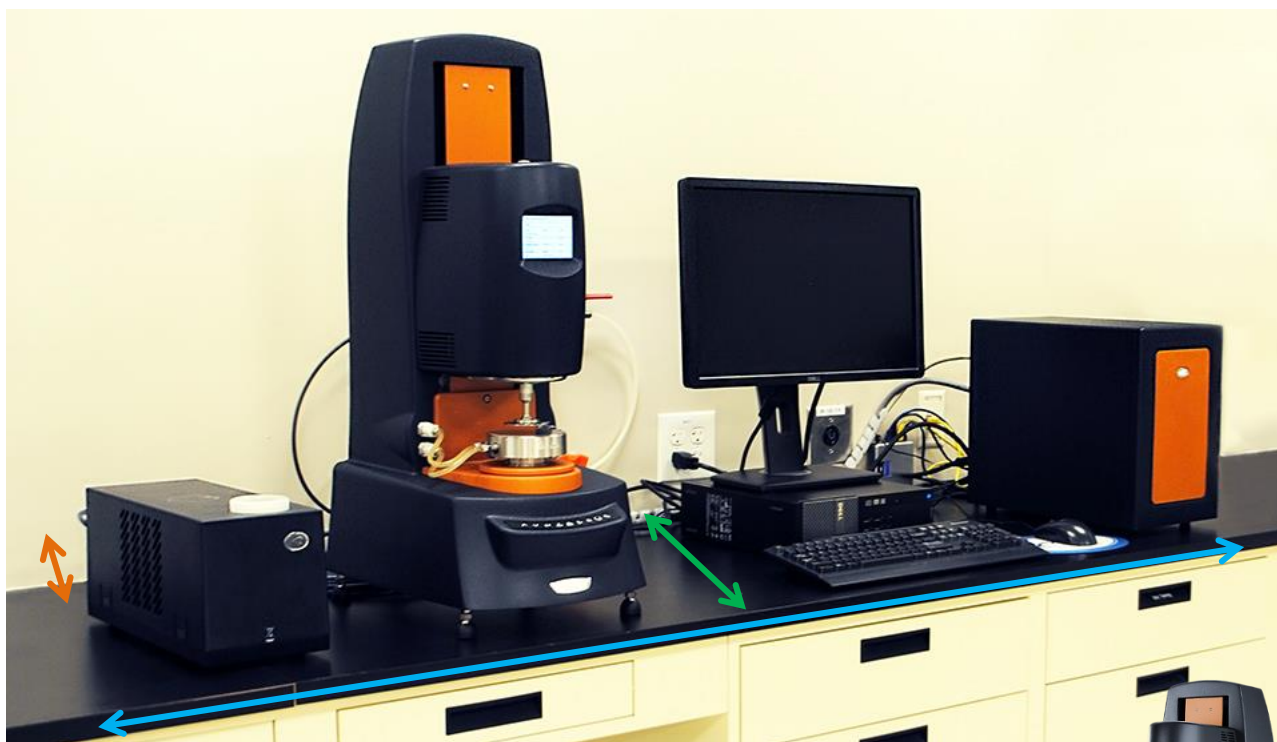


Bench length: 127 cm (50 in)

Marble table length: 60 cm (24 in)

Bench depth: 76 cm (30 in)

Marble table depth: 76 cm (30 in)


Distance from the wall: 30.5 cm (12 in) min.

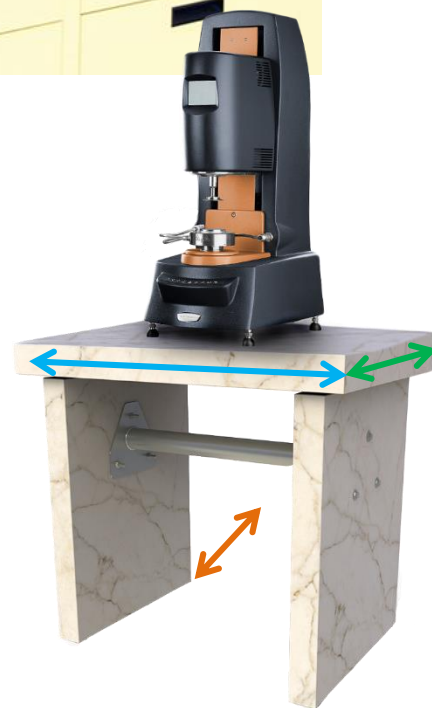

System Components



MAIN SYSTEM COMPONENTS

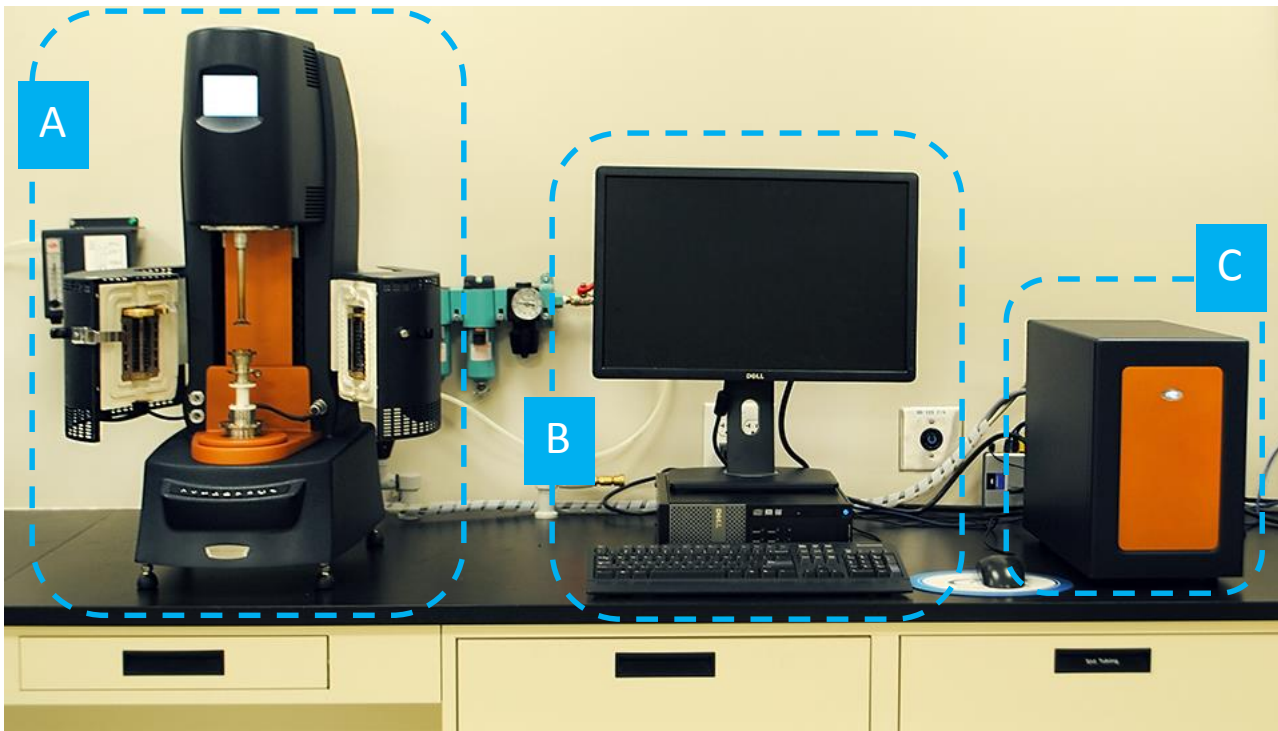

- A. Instrument
- B. Computer
- C. Electronics Control Module

Instrument Measurements



MAIN INSTRUMENT



Height: 76 cm (30 in)

Width: 32 cm (12.5 in)

Depth: 42 cm (16.5 in)

Weight: 32 kg (70.5 lbs)



ELECTRONICS CONTROL MODULE

Height: 48 cm (19 in)

Width: 26 cm (10 in)

Depth: 44 cm (17 in)

Weight: 14 kg (31 lbs)



Utility Requirements



POWER

110–230 VAC, 50/60 Hz, 1.4 kW



GAS

Air Bearing (Air or Nitrogen):

Gas Pressure	Compressed at 345–690 kPa gauge (50–100 psig)
Flow Rate	2 L/min
Dew point	-20°C or better
Conditions	<ul style="list-style-type: none">• Must be dry• Must be free from oil and dirt¹
Other	<ul style="list-style-type: none">• ¼ NPT female connection required for DHR main air supply (not provided)

¹Compressed Air Quality Requirements

Dew point	Ideal: -40°C Minimum: -20°C
Dirt particle	5µm
Oil including vapor	0.01 mg/m ³



WATER

Fluid circulator with cooling ability for Peltier and UHP temperature systems

Computer Requirements



HARDWARE REQUIREMENTS

Description	Requirement
Processor	<ul style="list-style-type: none">• Intel® Core™ 2 Duo or better• 2.93 GHz with 3 MB L2 cache
Memory	Required: ≥ 4 GB RAM Recommended: ≥ 8 GB RAM
Hard drive	≥ 80 GB free space <ul style="list-style-type: none">• 1 GB required for Full version of TRIOS• 600 MB required for Lite version of TRIOS (without Online help)
DVD	$\geq 48\times$ CD-ROM or DVD
Screen resolution	Required: 1280 x 1024 with 24-bit colors Recommended: 1920 x 1080 with 24-bit colors
Graphic memory	128 MB
Screen (LCD) size	Required: 19" or greater Recommended: 24" wide screen
USB II port	Required with ETC and Peltier Camera Viewer options, SALS accessory, and Automatic Asphalt calibration kit
Network card	Ethernet 10Base T/100 Base TX
Additional Ethernet card(s)	Necessary if connecting the instrument directly and access is needed to the Corporate LAN. Also required for Modular Microscope Accessory.
Ethernet Cabling	10/100BaseTX Ethernet hub/switch. Must be EIA-568B Category 5+ UTP
Client-Server Protocol	DHCP
Image Capture (Camera Option)	DirectX 9.0 or higher
Support for Custom Reporting	Rheology Advantage Navigator software only: Microsoft Word 97 or higher
Second Monitor	Recommended for SALS Accessory image viewing and Modular Microscope Accessory
TCP/IP ports used	<ul style="list-style-type: none">• TCP: 20010, 20011• UDP: 5050, 5056

Computer Requirements







SOFTWARE REQUIREMENTS

Item	TRIOS	Rheology Advantage
Operating System	<ul style="list-style-type: none">Windows 7, 8, 10 Ultimate, Enterprise & ProfessionalWindows Vista Business & UltimateHome version not supported	<ul style="list-style-type: none">Windows 7, 8 Ultimate, Enterprise & ProfessionalWindows Vista Business & UltimateHome version not supported
	Required: 32-bit or 64-bit version Recommended: 64-bit version	Required: 32-bit or 64-bit version
Browser	Internet Explorer	
Service Pack	Microsoft Operating System Service Pack	
Updates	Must be up-to-date	
Network	<i>TA Instruments is not responsible for resolving issues associated with connections to your corporate network.</i>	
Conflicts	<i>TA Instruments is not responsible for resolving hardware/software conflicts created by the addition of third party hardware or software to the computer.</i>	

Temperature Systems

The cooling rate and minimum temperature will depend on the source of cooling.

For accessory requirements not listed, refer to the DHR Accessory Requirements guide.

Accessory	Smart Swap Requirements	
Electrically Heated Plates (EHP)		<ul style="list-style-type: none"> Purge flow of 5 L/min (305 in³/min) inert gas Motor cooling gas flow of 10 L/min for temperatures above 250°C. Air pressure of 50–100 psig.
		<ul style="list-style-type: none"> Optional controlled cooling with <u>GCA</u>. Refer to the DHR Accessories Requirement guide for GCA requirements. Crash cooling pressure of 50–100 psig and a flow of ~2.5 scfm (70 L/min)
Environmental Test Chamber (ETC)		Purge gas flow rate should be 10 L/min (610 in ³ /min) at 206–690 kPa (30–100 psig)
Peltier Plate/ Peltier Concentric Cylinder		<ul style="list-style-type: none"> Recirculating water bath (not supplied) at 0.5 L/min (30.5 in³/min)



Circulator



Power



Cooling



Gas



LN₂



Fluid



Light



Hardware



Software



Temp














Lab



Customer

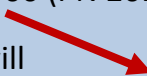Temperature Systems

Accessory	Smart Swap Requirements	
Upper Heated Plate (UHP)	 	<u>Option 1: Standard Cooling (temps above 10°C)</u> <ul style="list-style-type: none"> Circulation fluid: Requires continuous supply from mains tap or supplied fluid circulator Fluid cooling: Supply should be 5°C below the minimum required temperature at minimum required flow rate through the system of 0.5 L/min
	 	<u>Option 2: Standard Cooling Accessory</u> <ul style="list-style-type: none"> TA-supplied <u>Air Cooled Circulator</u> (PN 403209.901)  Recommended fluid: 100% distilled water with 1 oz. of TA conditioner added 
	 	<u>Option 3: Low Temperature Cooling Accessory–ThermoCube</u> <ul style="list-style-type: none"> TA-supplied <u>ThermoCube</u> model 10-300 (PN 201786.001)   Filling with 100% distilled water will damage the heat exchanger at 10°C and lower. Recommended fluid for UHP temp control set points <i>above</i> 10°C: 25% ethylene glycol, 75% distilled water with 1 oz. TA conditioner added Recommended fluid for UHP temp control set points <i>below</i> 10°C: 25% ethanol, 75% distilled water with 1 oz. TA conditioner added 



Circulator



Power



Cooling



Gas



LN₂



Fluid



Light



Hardware



Software



Temp







Lab



Customer

Temperature Systems

Accessory	Smart Swap Requirements	
Upper Heated Plate (UHP)		<u>Option 4: Low Temperature Cooling Accessory–Customer-supplied</u> <ul style="list-style-type: none"> Customer-supplied refrigerated and heating circulator and appropriate fluid (ie. silicone fluid)
		<div>  DO NOT USE WATER AS CIRCULATION FLUID </div> <ul style="list-style-type: none"> Supply: 5°C below the minimum required temperature at a minimum flow rate through the system of 0.5 L/min (12.2 in³/min)
		<u>Option 5: Low Temperature Cooling Accessory–Vortex</u> <ul style="list-style-type: none"> TA-supplied Vortex Cooler (PN 545809.901) Air: Clean, dry, oil-free, compressed air 200L/min at 552–690 kPa gauge (80–100 psig) Dew point: -30°C or better



Site Preparation Checklist



Discovery Hybrid Rheometer

	Sufficient bench space for instrument, computer, and Electronics Control Box <input type="checkbox"/> Length: 127 cm (50 in) <input type="checkbox"/> Depth: 76 cm (30 in)
	Instrument power is 110–230 VAC, 50/60 Hz, 1.4 kW
	Air Bearing Gas Pressure (air or nitrogen) <input type="checkbox"/> Pressure is 345–690 kPa (50–100 psig) <input type="checkbox"/> Flow rate is 2 L/min <input type="checkbox"/> Dew point is -20°C or better <input type="checkbox"/> ¼ NPT female connection to the main compressed air source
	Accessories used: <input type="checkbox"/> EHP <input type="checkbox"/> ETC <input type="checkbox"/> Peltier Place/Peltier Concentric Cylinder <input type="checkbox"/> Upper Heated Plate

I hereby acknowledge that all utility requirements have been met per the checklist above and that they will be ready at the agreed time of installation.

If all utility requirements are not met at the agreed time of installation, additional charges may be incurred for a return Service trip.

Customer *DD* / *MM* / *YYYY*

Company *City* *State* *Country*

Please send a signed copy of the completed checklist to your local Service representative.

TA Instruments Offices

For information on our latest products, contact information, and more, see our website at:
<http://www.tainstruments.com>.

To find your local TA Instruments office and contact information, visit
<http://www.tainstruments.com/contact/ta-directory/>

TA Instruments – Waters LLC
Corporate Headquarters
159 Lukens Drive
New Castle, DE 19720
USA

Telephone: 302-427-4000
Fax: 302-427-4001
Email: info@tainstruments.com